

WHEN YOU PUT PATIENTS' HEALTH FIRST, REMARKABLE THINGS FOLLOW

At Frankenman, we believe we're not so much in the business of making instruments but in the business of making surgeries more successful.

That's why we listen to doctors who've put our instruments to use. We read their reports. We invite them to think with us, to brainstorm, to help us develop new and innovative ideas.

That way we know we will be able to deliver the highest quality products to the medical profession both today and way into the future.



State-of-the-art manufacturing facilities.



CHEX
SURGICAL STAPLERS

INFORMED DESIGN.
SCRUPULOUSLY APPLIED.
METICULOUSLY PRODUCED.

OUR RANGE OF PRODUCTS INCLUDE

SINGLE USE CIRCULAR STAPLER FOR RECTAL PROLAPSE AND HAEMORRHOIDS

SINGLE USE RELOADABLE LINEAR STAPLER

SINGLE USE CURVED INTRALUMINAL CIRCULAR STAPLER

CS COMPACT™

SINGLE USE RELOADABLE LINEAR CUTTER STAPLER

SINGLE USE TRANSVERSE CUTTING LINEAR STAPLER



FOR MORE INFORMATION OR TO VIEW OUR ENTIRE PRODUCT RANGE PLEASE VISIT OUR WEBSITE AT WWW.FRANKENMAN.HK



FRANKENMAN INTERNATIONAL LTD
SUITE B 13/F WING TAT COMMERCIAL BUILDING
121-125 WING LOK STREET SHEUNG WAN HONG KONG
T +852 3106 3035 F +852 3585 0519
E info@frankenman.hk

www.frankenman.hk
Chex and Chex Healthcare are registered trademarks of Frankenman International Limited



"UNRIVALED VISUALIZATION -
A NEW PERSPECTIVE ON THE PROCEDURE"

INFORMED DESIGN.
SCRUPULOUSLY APPLIED.
METICULOUSLY PRODUCED.



HIGH VOLUME STAPLED
HAEMORRHOIDOPEXY
with the CPH34 HV

(PATENT PENDING)

CHEX
SURGICAL STAPLERS

CPH34HV - THE FIRST AND ONLY DEVICE DESIGNED FOR STAPLED HAEMORRHOIDAL PROLAPSECTOMY

FEATURES AND BENEFITS

- The Chex™ CPH34HV remarkable transparent housing gives the surgeon a new perspective on the procedure that was not possible until now
- Lowered tissue stop to facilitate 'pull through' of tissue into housing
- The largest volume staple housing on the market today allowing the resection of more tissue when clinically required



In surgical treatment of haemorrhoidal disease, and especially when combined with ODS, a key surgical objective is the correction of internal rectal prolapse.

Until the introduction of the Chex™ CPH range from Frankenman International Ltd staplers used in the treatment of haemorrhoids had evolved from conventional surgical staplers, designed for conventional surgical anastomosis making their use in the treatment of associated haemorrhoidal prolapse less than ideal.

The Chex™ CPH34HV is the first device specifically designed to maximise the volume of the staple housing whilst maintaining the essential design characteristics of, what has been proven scientifically and in clinical use, to be the most effective device in use today for Stapled haemorrhoidopexy.


The Chex™ CPH34HV a remarkable step forward in the treatment of haemorrhoids, ODS and associated Rectal Prolapse.



SIDE VIEW – ACTUAL SIZE

KEY FEATURES & SPECIFICATIONS

- 1 STAINLESS STEEL ANVIL SHAFT TO MAXIMIZE STABILITY DURING FIRING AND CLAMPING
- 2 HIGH VOLUME HOUSING
- 3 TRANSPARENT STAPLE HOUSING AND STAPLE PUSHERS GIVING MAXIMUM VISUAL FEEDBACK
- 4 SUTURE CONDUITS ALLOWING BETTER ACCESS AND ADDITIONAL OPTIONS IN TISSUE MANIPULATION

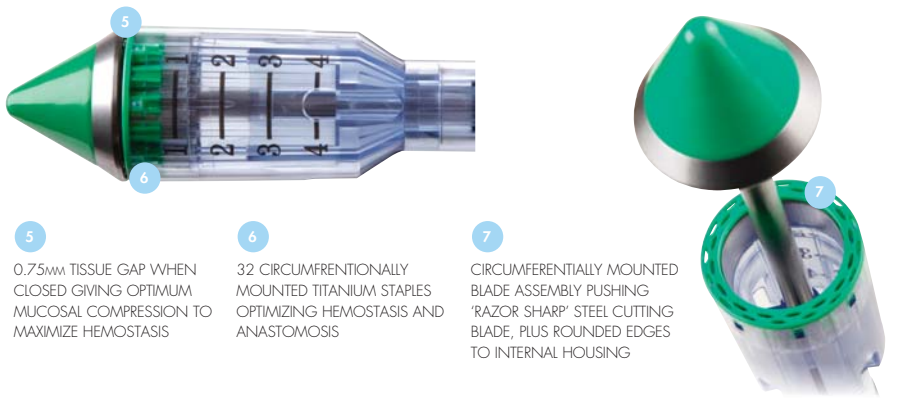


- 5 0.75mm TISSUE GAP WHEN CLOSED GIVING OPTIMUM MUCOSAL COMPRESSION TO MAXIMIZE HEMOSTASIS
- 6 32 CIRCUMFERENTIALLY MOUNTED TITANIUM STAPLES OPTIMIZING HEMOSTASIS AND ANASTOMOSIS
- 7 CIRCUMFERENTIALLY MOUNTED BLADE ASSEMBLY PUSHING 'RAZOR SHARP' STEEL CUTTING BLADE, PLUS ROUNDED EDGES TO INTERNAL HOUSING

PRODUCT CODE	CARTRIDGE COLOUR	STAPLE ROWS	OUTSIDE DIAMETER	CUTTING DIAMETER	STAPLE LEG LENGTH	STAPLE CROWN LENGTH	STAPLE COUNT	CLOSED STAPLE HEIGHT
CPH34HV	GREEN WITHIN TRANSPARENT HOUSING	2	34.5mm	26mm	4.2mm	3.8mm	32	0.75mm – 1.5mm

HIGHEST VOLUME HOUSING

The highest Volume of tissue resected in comparison market leading products when measure in Weight (g), Volume (ml) and cubic millimetres.

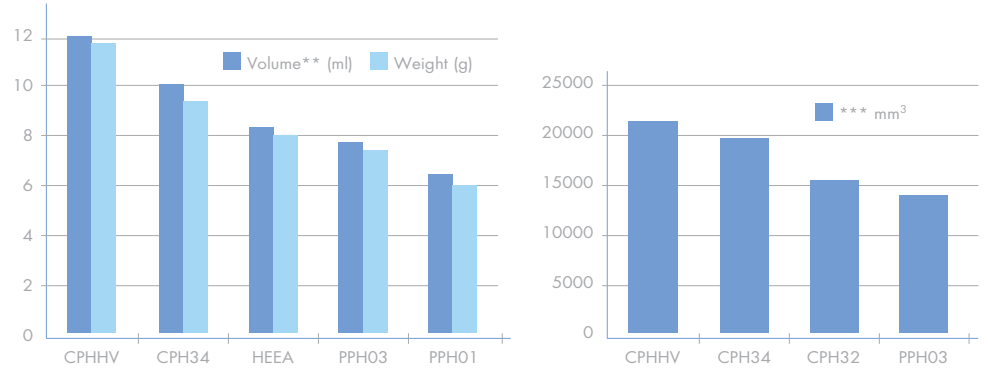

*Technological Improvements in the Treatment of Hemorrhoids and Obstructed Defecation Syndrome. *Experimental Study, Giuliano Reboa*, Marco Gipponi*, Tommaso Testa* and Francesca Lantieri***; *U.O.C. General Surgery, Colo-Rectal Unit - "San Martino" Hospital, Genoa - Italy; **Health Science Dept, Biostatistic Unit - University of Genoa, School of Medicine. In Vivo January 2011 25:129-135

SET CONTAINS: CIRCULAR STAPLER FOR PROLAPSE AND HAEMORRHOIDS(CPH), SUTURE THREADER, ANAL DILATORS (STANDARD AND "BUTTERFLY"), ANAL DILATOR INNER LINING AND PURSE STRING SUTURE ANOSCOPE